



## USING THE QUESTIONNAIRES IN THE WIDEX BABY440 FITTING WITH COMPASS V5.1

This bulletin describes the questionnaires in the WIDEX BABY440 fitting using Compass V5.1. The description refers to the *Parent questionnaire* and *Hearing checklist* in the Client section and Finalise section. Please refer to Compass in order to see details of the screen. Remember that Compass offers several help functions to explain the different options in the program: Tooltips, Solution guide, Using this window panels and the Help Manual.

With ChildFit, two new questionnaires have been introduced in Compass V5.1. – The *Parent questionnaire* and the *Hearing checklist*. The questionnaires are used to validate the fitting of the hearing aid and to follow the auditory, language and speech development of the child. You find the questionnaires in the Client window and Finalise window in Compass.

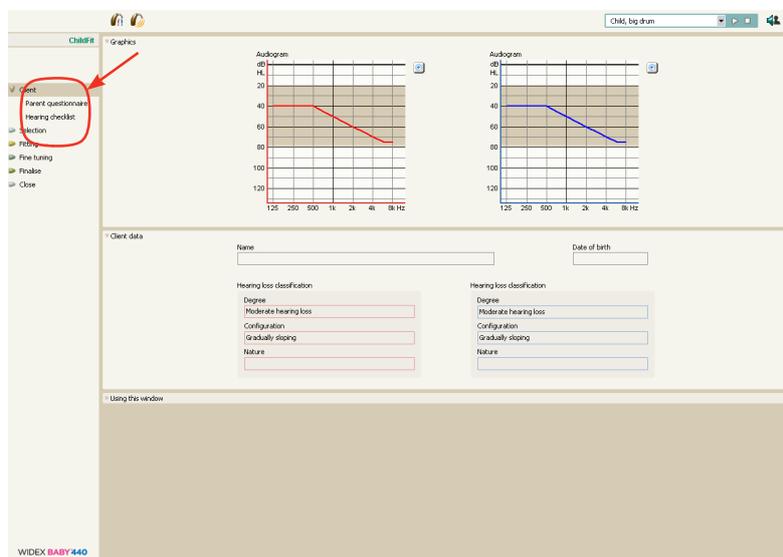

Finding the questionnaires under the Client window

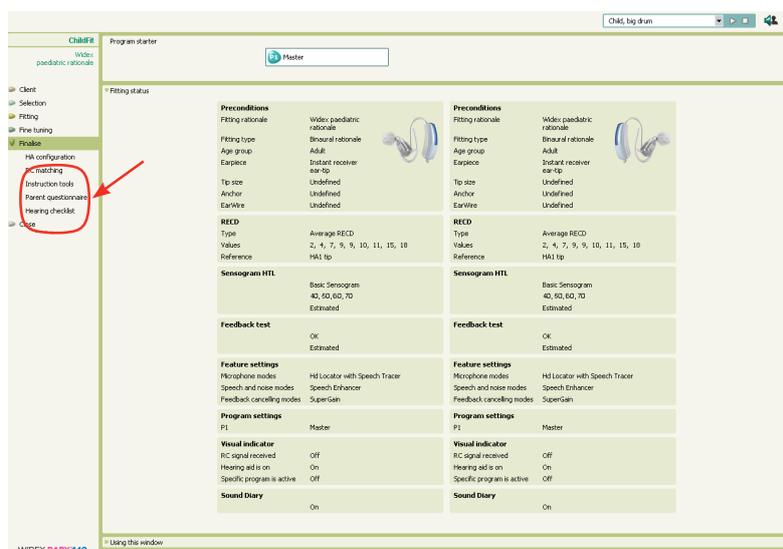

Finding the questionnaires under the Finalise window

## USING THE PARENT QUESTIONNAIRE

The *Parent Questionnaire* is used to validate the fitting of the hearing aid and contains 16 questions divided into four parts. This questionnaire contains questions concerning the physical fit, daily use and reactions during the first period the child wears hearing aids.

The screenshot shows the 'Parent questionnaire' window in the Compass v. 5.1 software. At the top, there is a client data section with fields for Name, Date of birth, and Age group. Below this, the questionnaire is displayed with a 'Select questions' dropdown menu set to 'All (1-16)'. The questionnaire consists of four questions:

1. How many hours are the hearing aids worn each day?  
 No info  
 All day  
 4-8 hours  
 1-4 hours  
 Less than 1 hour
2. Do the hearing aids whistle when worn?  
 No info  
 Never  
 Sometimes  
 Often  
 Always
3. Does the earpiece fit well in your child's right ear?  
 No info  
 Yes  
 No  
If no, please specify the reason:  
 Slips out of ear  
 Loose fit  
 Irritation in the ear  
 Allergic reaction  
 Other
4. Does the earpiece fit well in your child's left ear?  
 No info  
 Yes  
 No  
If no, please specify the reason:  
 Slips out of ear  
 Loose fit  
 Irritation in the ear  
 Allergic reaction  
 Other

The Parent questionnaire in Compass v. 5.1.

When first opening the *Parent questionnaire*, all 16 questions are displayed. Use the drop-down list above the questionnaire to choose which part you want displayed.

Some of the questions are followed by radio buttons where only one answer can be entered. Others are followed by a number of tick boxes allowing more than one answer, and finally some questions are followed by a text field for comments to be entered.

You can use the *Parent questionnaire* in the following way:

1. Open your client record from your database
2. Open the *Parent questionnaire* window
3. Check that the information at the top of the window is correct.
4. Answer as many of the questions as you can.
5. If desired, enter comments for the questions that are equipped with a comment field

## USING THE HEARING CHECKLIST

In the *Hearing checklist*, information about the child and the child's development is entered. This can help the dispenser follow the individual child's auditory, language and speech development, and thereby decide whether adjustments of the fitting are needed.

The hearing checklist in Compass v. 5.1.

The *Hearing checklist* contains a list of auditory, language and speech milestones categorised into six levels. Each level corresponds to the development of milestones seen in normal hearing children, in the following way:

Level	Age
Level 1	0-3 months
Level 2	3-6 months
Level 3	6-12 months
Level 4	12-18 months
Level 5	18-24 months
Level 6	24-36 months

You should however expect that the hearing impaired child's auditory development differs from that of normal-hearing children and therefore the normal age range should only be considered a guideline.

When first opening the *Hearing checklist*, the first level is displayed. Use the drop-down list above the checklist to change to another level. If *All* is selected from the drop-down list, it is possible to scroll through all the milestones.

For each of the milestones information can be entered about when the child acquired the skill, and you can enter whether the child uses the skill *Always*, *Sometimes* or *Never* to a given test. It is also possible to indicate the date on which the child acquired a certain skill by using the calendar icon .

You can use the *Hearing checklist* in the following way:

1. Open your client record from your database.
2. Open the *Hearing checklist* window
3. Check that the information at the top of the window is correct.
4. Select the level that corresponds to your client's age.
5. Answer as many of the items in the checklist as you can. For each item, select the appropriate answer. You can only enter one answer for each item in the list.

You may need to go back to a previous level to enter more information. Use the drop-down list above the checklist to move to another level.

## STORING OF INFORMATION AND FOLLOW UP

Information entered into the questionnaires on the screen is saved in the client files in the client database. This is helpful for the audiologist during follow-up visits in order to see any progress in the fitting of the hearing aids or in the development of the child. At the top of the questionnaires you can see when the questionnaire was first created and when it was updated.

## PRINTING THE QUESTIONNAIRES

If you have Microsoft Word installed on your computer, you can also print the questionnaire. Before using the Print icon in the window, choose what you want to print. You can choose between printing an empty questionnaire for the parents to fill out at home, or you can print a questionnaire containing the answers that you have already filled out in Compass.