

WIDEX MRR4D CHARGER

**MAIN FEATURES**

- Compact and portable.
- Can be powered from multiple power sources.
- Intuitive use supported by light indicators displaying charge level of the hearing aid.

DESCRIPTION

The WIDEX sRIC charger is designed to charge two WIDEX MRR4D hearing aids.

The charger has a charging well for each hearing aid, and two light indicators that show the charging status of the individual hearing aid. Each hearing aid can be charged individually, depending on the state of charge of the hearing aid. This means that if one hearing aid is fully charged and placed in the charger, it will not charge. However, if the other hearing aid is not charged, it will charge.

PROPERTIES	
Colour:	Black and warm platinum (metallic)
Dimensions:	Height: 51mm Width: 77mm Depth: 27 mm
Weight:	approx. 85 g
Capacity:	Holds a set of hearing aids. Can charge the hearing aids individually. Charging time and indications see Instruction for Use for the WIDEX MRR4D Charger
Operating temperature, humidity and atmospheric pressure ranges:	0°C (32°F) to 30°C (86°F), 10%-95% rH and 750 to 1060 mbar.
Storage and transportation:	-20°C (-4°F) to 50°C (122°F), 10%-95% rH and 750 to 1060 mbar

POWER SUPPLY	
Input	Type: DC, stabilized 5V +/- 5% Min current: 500 mA Connector: micro USB receptacle Different power supplies can be used (see User Instructions). The hearing aid can become warm during charging.


DISPOSAL

- Do not dispose of the charger (new or used) along with ordinary household waste.
- The charger must be disposed of as electronic or chemical waste, according to local regulations.

The Charger at a glance

The charger has a charging well for each hearing aid, and two light indicators that show the charging status of the hearing aids.

1. Micro USB port
2. Charging wells
3. Light indicators

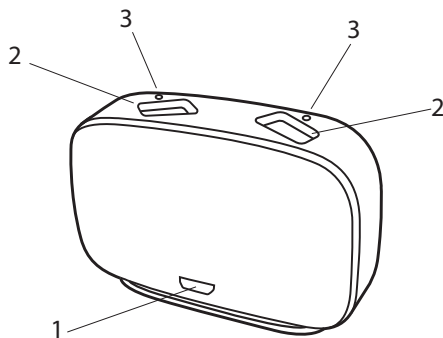

IN THE BOX

In the presentation box you will find:

- One charger
- One power adaptor
- One USB cable
- User instructions